

Lifeline




MedReady[®] Medication Dispenser

A simple, affordable and automated way to manage multiple medications with 24/7 monitoring

Medication errors are common and costly


You rely on your medications to keep you healthy. But complicated medication schedules can lead to mistakes: missing doses, taking incorrect amounts, or taking pills at the wrong time. These mistakes could lead to unnecessary doctor or hospital visits.

The MedReady Medication Dispenser* helps you to avoid medication mistakes so you can live independently in the comfort of your own home. Our automated medication dispenser with 24-hour monitoring is designed to help you to take the right dose at the prescribed time.

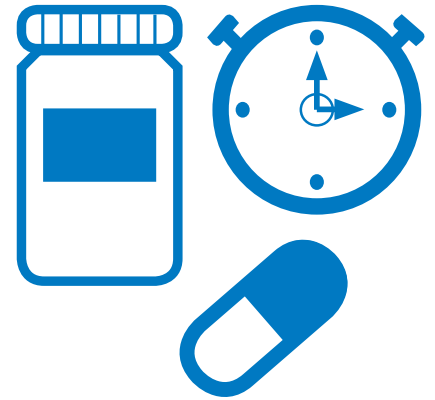

How it works

- 1. Fill:** Using the dosing templates (1 to 4 doses a day) add medications from dose 1 then confirm when completed.
- 2. Take medication:** The pill tray rotates to the scheduled dose and the audible alarm sounds. Open the medication compartment and retrieve the pills.
- 3. Alerts:** Caregiver receives an SMS text or email¹ if medication has not been taken or when the dispenser needs troubleshooting.
- 4. Lifeline Response Centre support:** Our specially trained Associates will assist the caregiver and can send the right help if needed.

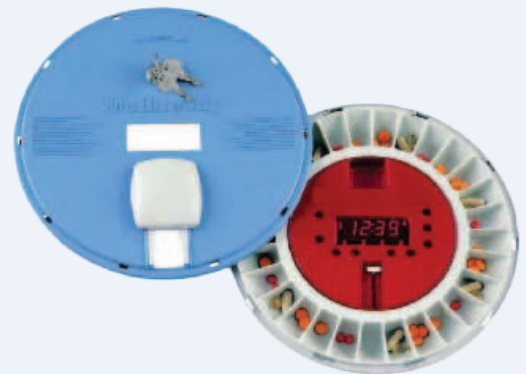

Is it time for a better way to manage your medications?

Use this self-assessment as a guide to determine if a medication management device is right for you.

Questions (Check the box if your answer is “Yes”)


1. Are you taking three or more medications each day?
2. Do you take medicine two or more times per day?
3. In the past month, have you missed medications because you forgot to take them?
4. Do you have a chronic condition or disease, such as: heart disease, COPD, osteoporosis, diabetes, or cancer?
5. Do you require assistance with daily activities, such as taking medication, bathing, dressing, or cooking?
6. Would a highly reliable service for dispensing medications accurately and on time provide peace of mind for you or your loved ones?
7. Is it important for you to continue living independently in your current home?
8. Do you have concerns about complications that could lead to hospitalizations if you don't take your medicine correctly?
9. Are your daily activities restricted by your medication regimen?


See the chart at the right to review the assessment for your score.


6 – 9 Yes answers:	This high score indicates that you may be at serious risk for medication mistakes and resulting complications.
3 – 5 Yes answers:	You could be at some risk for medication mistakes and resulting complications.
1 – 2 Yes answers:	Ask yourself two important questions: <ol style="list-style-type: none">1. If I don't take my medicine correctly, am I at risk for further complications?2. If I don't take my medicine correctly, am I at risk of having to leave my home so that someone can help me take them?


If the answer is “yes” to any of these questions, it might be a good idea to discuss an automated medication management service with your family and healthcare provider.

Features and Benefits

 **Easy to use** - Dispenses pre-loaded doses at the scheduled time². Simply press open the compartment to receive the medication.

 **Monitored by Lifeline's Response Centre** - Caregiver is alerted when a dose is missed by SMS text or email. The Response Centre is also available to support the caregiver when needed.

 **Highly reliable** - Visual and audible alarms indicate when it is time to take the medication and allows for greater adherence.

 **Back up battery** - Designed to distribute every dose on time, even in a power outage, thanks to our reliable backup battery³ protection.

Why choose Lifeline?

- ▶ Lifeline is the #1 medical alert service in Canada⁴
- ▶ Caring, friendly Lifeline Response Associates are ready to help 24/7
- ▶ Canadian round-the-clock monitoring
- ▶ No long-term contract
- ▶ Our closed-loop protocol means we follow-up to ensure help arrives
- ▶ Proven fall detection technology available

MedReady® is a registered trademark of MedReady, Inc.

* This service is only available at participating Lifeline programs.

1. The subscriber and/or caregiver require a smart phone in order to enroll in this service.

2. Certain limitations subject to third party cellular provider availability and coverage. Signal range may vary.

3. Back up power battery can last up to 72 hours.

4. Based on number of subscribers.

1-800-LIFELINE
1-800-543-3546